



366PR-05

G-Biosciences ♦ 1-800-628-7730 ♦ 1-314-991-6034 ♦ technical@GBiosciences.com

A Geno Technology, Inc. (USA) brand name

PAGEmark™ Blue PLUS Protein Markers

Prestained Protein Ladder for
Molecular Weight Determination by Electrophoresis

(Cat. # 786-418)



think proteins! think G-Biosciences www.GBiosciences.com

INTRODUCTION

A blue protein standard with 12 prestained proteins covering a wide range of molecular weights:


- 10-240kDa in Tris-Glycine buffer
- 9-235kDa in Bis-Tris (MOPS) buffer
- 9-235kDa Bis-Tris (MES) buffer

Proteins are covalently coupled with a blue chromophore and two reference bands, 25kDa and 72kDa, are enhanced in intensity for easier reference.

The PAGEmark™ line do not require boiling or denaturing, simply load on gel, and are readily transferred to Western blot membranes (PVDF and nitrocellulose).

ITEM(S) SUPPLIED

Cat. #	Product	Loads
786-418	PAGEmark™ Blue PLUS	100 Loads with 5µl/well

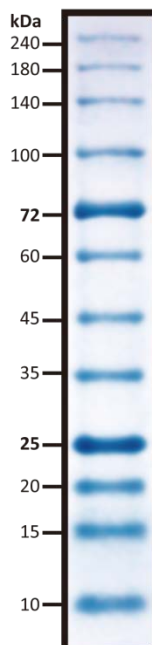

4-20% Tris-Glycine

STORAGE CONDITIONS

All PAGEmark™ markers are shipped at ambient temperature.

Upon arrival, store all PAGEmark™ markers at -20°C. Markers are stable for up to a month at 4°C.

IMPORTANT INFORMATION

- We highly recommend that the larger volumes markers are aliquoted into smaller fractions. This will help prevent contamination and repeated freeze-thawing may damage the product.
- The markers are not suitable for native polyacrylamide gel electrophoresis
- For silver staining, dilute the markers approximately 1:10 in reducing sample loading buffer.
- ***These markers are ready-to-use and do not require boiling, dilution or addition of reducing agents before loading***

PROTOCOL

1. If a precipitate is visible, warm to 40°C before use.
2. Loading volumes:
 - a. 3-5µl/well on mini gels for clear visualization during electrophoresis.
 - b. 1.2-2.5µl per well for general Western transfers
 - c. Apply more for thicker (>15mm) or larger gels

MOLECULAR WEIGHT ESTIMATION GUIDE

Color	Tris-Glycine	Bis-Tris (MOPS)	Bis-Tris (MES)
Blue	240	235	235
	180	170	170
	140	130	130
	100	93	93
	72	68	70
	60	53	53
	45	41	42
	35	30	30
	25	22	23
	20	18	18
	15	14	14
	10	9	9

kDa

240

180

140

100

72

60

45

35

25

20


15

10

4-20% Tris-Glycine

RELATED PRODUCTS

Download our Protein Electrophoresis and Western Blotting Handbooks.

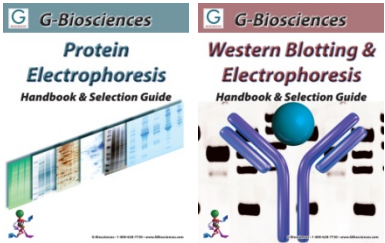

<http://info.gbiosciences.com/complete-protein-electrophoresis-handbook/>

<http://info.gbiosciences.com/complete-western-blot-handbook--selection-guide/>

For other related products, visit our website at www.GBiosciences.com or contact us.

Last saved: 12/17/2015 CMH



www.GBiosciences.com