



366PR-05

G-Biosciences ♦ 1-800-628-7730 ♦ 1-314-991-6034 ♦ technical@GBiosciences.com

A Geno Technology, Inc. (USA) brand name

PAGEmark™ Tricolor PLUS Protein Markers

Prestained Protein Ladder for
Molecular Weight Determination by Electrophoresis

(Cat. # 786-419)



think proteins! think G-Biosciences www.GBiosciences.com

INTRODUCTION

A ready-to-use three-color protein standard with 13 prestained proteins covering a wide range of molecular weights:


- 5-245kDa in Tris-Glycine buffer
- 3.5-235kDa in Bis-Tris (MOPS) buffer
- 3.5-235kDa Bis-Tris (MES) buffer).

Proteins are covalently coupled with a blue chromophore except for two reference bands. A green chromophore at 25kDa and a red chromophore at 75kDa.

The PAGEmark™ line do not require boiling or denaturing, simply load on gel, and are readily transferred to Western blot membranes (PVDF and nitrocellulose).

ITEM(S) SUPPLIED

Cat. #	Product	Loads
786-419	PAGEmark™ Tricolor PLUS	100 Loads with 5µl/well

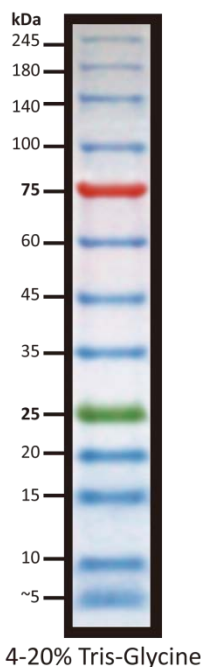

STORAGE CONDITIONS

All PAGEmark™ markers are shipped at ambient temperature. Upon arrival, store all PAGEmark™ markers at -20°C. Markers are stable for up to a month at 4°C.

IMPORTANT INFORMATION

- We highly recommend that the larger volumes markers are aliquoted into smaller fractions. This will help prevent contamination and repeated freeze-thawing may damage the product.
- The markers are not suitable for native polyacrylamide gel electrophoresis
- For silver staining, dilute the markers approximately 1:10 in reducing sample loading buffer.
- ***These markers are ready-to-use and do not require boiling, dilution or addition of reducing agents before loading***

PROTOCOL

1. If a precipitate is visible, warm to 40°C before use.
2. Loading volumes:
 - a. 3-5µl/well on mini gels for clear visualization during electrophoresis.
 - b. 1.2-2.5µl per well for general Western transfers
 - c. Apply more for thicker (>15mm) or larger gels


MOLECULAR WEIGHT ESTIMATION GUIDE

Color	Tris-Glycine	Bis-Tris (MOPS)	Bis-Tris (MES)
Blue	245	235	235
Blue	180	170	170
Blue	140	130	130
Blue	100	93	93
Red	75	70	72
Blue	60	53	53
Blue	45	41	42
Blue	35	30	30
Green	25	22	23
Blue	20	18	18
Blue	15	14	14
Blue	10	9	10
Blue	~5	3.5	3.5

4-20% Tris-Glycine

RELATED PRODUCTS

Download our Protein Electrophoresis and Western Blotting Handbooks.

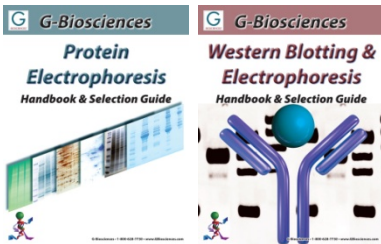

<http://info.gbiosciences.com/complete-protein-electrophoresis-handbook/>

<http://info.gbiosciences.com/complete-western-blot-handbook--selection-guide/>

For other related products, visit our website at www.GBiosciences.com or contact us.

Last saved: 12/17/2015 CMH



www.GBiosciences.com