



PRODUCT INFORMATION

Prestained Protein Molecular Weight Marker

Pub. No. MAN0011770

Rev. Date 19 June 2019 (Rev. C.00)

Components	#26612
Prestained Protein Molecular Weight Marker	2 x 250 µL

Store at -20 °C

www.thermofisher.com

www.thermoscientific.com/pierce

For Research Use Only. Not for use in diagnostic procedures.

Introduction

The Prestained Protein Molecular Weight Markers are a mixture of proteins ranging from 20 kDa to 120 kDa producing six well-defined blue bands by SDS-PAGE. The protein marker mix is conveniently packaged and ready to use with no heating, diluting or additional reducing agent necessary. Lot-to-lot variation of the apparent molecular weight of prestained proteins is ~5 %.

Storage Buffer: 62.5 mM Tris•H₃PO₄ (pH 7.5 at 25 °C), 1 mM EDTA, 2 % (w/v) SDS, 10mM DTT, 1.5 mM NaN₃, and 33 % (v/v) glycerol.

Important Product Information


- Do not boil the marker mix.
- Do not use the marker to determine the native molecular weight of proteins in native polyacrylamide gel electrophoresis.
- The mobility of prestained proteins can vary in different SDS-PAGE buffer systems; however, they are suitable for approximate molecular weight determination.
- For precise protein MW determination use the Thermo Scientific™ PageRuler™ Broad Range Unstained Protein Ladder (#26630).

Procedure for Use in Polyacrylamide

Gel Electrophoresis

1. Thaw the marker at room temperature. Do not boil the protein marker.
2. Mix gently and thoroughly to ensure that the solution is homogeneous.
3. Load an appropriate volume of the marker mix onto the gel.
 - Mini-gel: 5 μ L per well (0.75-1.0 mm thick) or 10 μ L per well (1.5 mm thick)
 - Large gel: 10 μ L per well (0.75-1.0 mm thick) or 20 μ L per well (1.5 mm thick)
4. Store the unused protein marker at -20 °C.

Prestained Protein Molecular Weight Marker

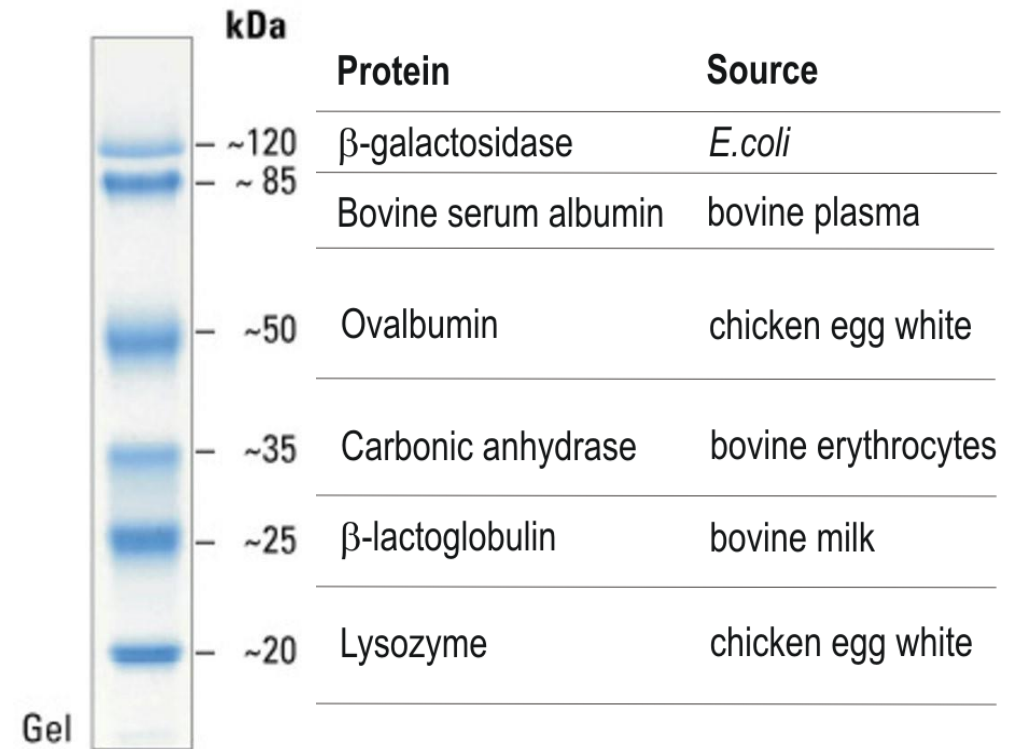

8-16% Tris-glycine SDS-PAGE

General References

Burnette, W.N. (1981). "Western blotting": electrophoretic transfer of proteins from sodium dodecyl sulfate – polyacrylamide gels to unmodified nitrocellulose and radiographic detection with antibody and radioiodinated protein A. *Anal Biochem* 112(2):195-203.

Kurien, B.T. and Scofield, R.H. (2003). Protein blotting: a review. *J Imm Meth* 274:1-15.

Laemmli, U.K. (1970). Cleavage of structural proteins during the assembly of the head of bacteriophage T4. *Nature* 227:680-5.

Towbin, H., et al. (1979). Electrophoretic transfer of proteins from polyacrylamide gels to nitrocellulose sheets: procedure and some applications. *Proc Natl Acad Sci USA* 76:4350-4.

PRODUCT USE LIMITATION

This product is developed, designed and sold exclusively *for research purposes and in vitro use only*. The product was not tested for use in diagnostics or for drug development, nor is it suitable for administration to humans or animals.

Please refer to www.thermofisher.com for Material Safety Data Sheet of the product.

© 2019 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries.