

Blue Plus[®] II Protein Marker (14-120 kDa)

Please read the manual carefully before use

Cat. No. DM111

Storage: at -18°C or below for two years

Concentration

about 2 µg/5 µl for each band

Description

Blue Plus[®] II Protein Marker is composed of eight prestained proteins ranging from 14 to 120 kDa. The proteins of 50 kDa and 120 kDa bands are covalently coupled to orange dye. The protein of 14 kDa band is covalently coupled to yellow dye. The other five bands are covalently coupled to blue dye. This prestained protein marker is designed for monitoring the progress of SDS-polyacrylamide gel electrophoresis, for assessing transfer efficiency onto PVDF or NC membranes. Clear colored bands can be visible on PVDF or NC membranes. The orange and bright yellow bands make it easy to determine the direction of the transfer.

Storage buffer

100 mM Tris-HCl (pH 6.8), 5 mM EDTA, 10 mM DTT, 10% Glycerol, 1% SDS, 0.01% phenol red.

Protocol

- Ready-to-use, direct load on gels without heating and reducing agents.
- Thoroughly thaw and mix it. Use 5 µl/well for 1 mm thick mini-gel, 10 µl/well for larger gel.


Experimental result

After 12% Tris-Glycine SDS gel electrophoresis, transfer to PVDF membrane.

Notes

The **size** of the protein covalently coupled to dye will vary during **gel electrophoresis at different concentrations**, and it is only used as a reference when judging the molecular weight of the target protein.

Blue Plus[®] II Protein Marker
(14-120 kDa)

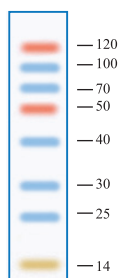

12% Tris-glycine SDS gel (5 µl/well)

BioRad Mini Electrophoresis System, 200 V 50 minutes

For research use only, not for clinical diagnosis.

Service telephone +86-10-57815020

Service email custserv@transgenbiotech.com

