


BLUeye Prestained Protein Ladder (A318979)


Specifications:

Name:	BLUeye Prestained Protein Ladder
Description:	Pre-stained protein ladder. Covers the size range 10-245 kDa (Tris-glycine-SDS running buffer). Number of bands = 12. Green band at approximately 25 kDa. Red band at approximately 75 kDa. Designed for monitoring SDS-PAGE separation and WB transfer efficiency. Thaw at room temperature or 37-40°C for a few minutes prior to use to dissolve precipitated solids. Do not boil. Do not heat, dilute, or add reducing agent before loading. Mix thoroughly.
Applications:	SDS-PAGE, WB
Formulation:	20 mM Tris-Phosphate pH7.5, 2% SDS, 0.2 mM Dithiothreitol, 3.6 M Urea and 15% (v/v) Glycerol.
Product Form:	Liquid
Storage:	Shipped at 4°C. Store at 25°C for up to 2 weeks, 4°C for up to 3 months, and -20°C for long-term storage.
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.

Images:

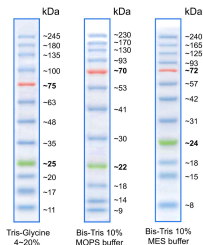

The prestained protein ladder contains 12 proteins in the range 10-245 kDa. Green and red reference bands are included at 25 and 75 kDa, respectively.

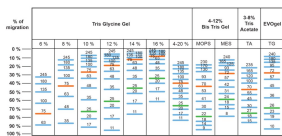

The prestained protein ladder contains 12 proteins in the range 10-245 kDa. Green and red reference bands are included at 25 and 75 kDa, respectively. The apparent molecular weight of each protein was determined by calibration against an unstained protein ladder for each electrophoresis condition.