


PiNK Plus Prestained Protein Ladder (A318977)

Specifications:

Name:	PiNK Plus Prestained Protein Ladder
Description:	Pre-stained protein ladder. Covers the size range 10-175 kDa (Tris-glycine-SDS running buffer). Number of bands = 11. Blue stained bands at approximately 10, 40 and 90 kDa. Designed for monitoring SDS-PAGE separation and WB transfer efficiency. Thaw at room temperature or 37-40°C for a few minutes prior to use to dissolve precipitated solids. Do not boil. Do not heat, dilute, or add reducing agent before loading. Mix thoroughly.
Applications:	SDS-PAGE, WB
Formulation:	20 mM Tris-Phosphate pH7.5, 2% SDS, 1 mM Dithiothreitol, 4.8 M Urea and 12% (v/v) Glycerol.
Product Form:	Liquid
Storage:	Shipped at 4°C. Store at 25°C for up to 2 weeks, 4°C for up to 3 months, and -20°C for long-term storage.
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.

Images:

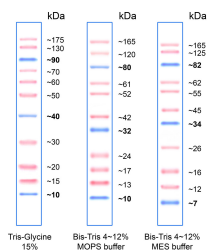

The prestained protein ladder contains 11 proteins that resolve into sharp bands in the range 10-175 kDa.



The prestained protein ladder contains 11 proteins that resolve into sharp bands in the range 10-175 kDa. The apparent molecular weight of each protein was determined by calibration against an unstained protein ladder for each electrophoresis condition.