

ProteinRuler® II (12-120 kDa)

Please read the manual carefully before use

Cat. No. DR201

Version No. Version 1.0

Storage: at -18°C or below for two years

Concentration

2 µg/5 µl for 12 kDa band; 1 µg/5 µl for 50 kDa band; 0.5 µg/5 µl for each of other bands

Description

ProteinRuler® II is composed of eight unstained recombinant proteins (20 kDa, 30 kDa, 40 kDa, 50 kDa, 60 kDa, 80 kDa, 100 kDa, 120 kDa) and one blue prestained recombinant protein (12 kDa). The prestained band allows monitoring electrophoresis. The 50 kDa band has doubled intensity to serve as a reference band. After SDS-PAGE electrophoresis, nine distinct bands can be visualized with Coomassie Brilliant Blue staining.

- MW range from 12 to 120 kDa.

Storage buffer

50 mM Tris-HCl (pH 6.8), 5 mM EDTA, 10 mM DTT, 10% glycerol, 1% SDS, 0.01% phenol red


Protocol

- Ready-to-use, direct load on gels without heating and reducing agents.
- Thoroughly thaw and mix it. Use 5 µl/well for 1 mm thick mini-gel, 10 µl/well for larger gel.

Electrophoresis condition

For BioRad Mini Electrophoresis installation, electrophoresis at 200 V, 50 minutes.

ProteinRuler® II
(12 kDa-120 kDa)

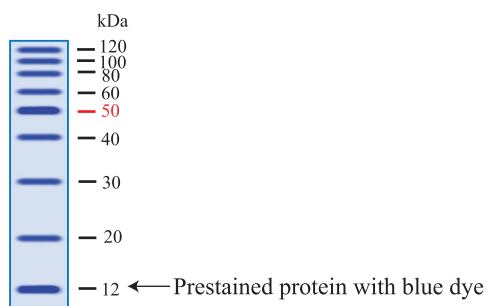

5 µl of *ProteinRuler*® II on 12% Tris-glycine SDS gel, stained with Coomassie Brilliant Blue.

For research use only, not for clinical diagnosis.

Service telephone +86-10-57815020

Service email custserv@transgenbiotech.com

